



ETS 5532/5533 Environmental Chamber

Operating Manual

Document D01042 Revision D

This manual covers the ETS 5532 chamber and the ETS 5533 chamber. Operating systems and set-up are identical for these two chambers. The essential difference is the physical depth of the environmental chamber – 6” larger in the case of 5533.

Table of Contents

I.	Important Safety Information	3
II.	Description of Contents	4
III.	Setup Guide	5
IV.	Quick Start Guide	9
V.	Operation Guide	10
VI.	Calibration and Maintenance	14
VII.	Troubleshooting	16
VIII.	System Specifications	18
IX.	Warranty	20

I. Important Safety Information

SAFETY INSTRUCTIONS

The equipment described in this Manual is designed and manufactured to operate within defined design limits. Any misuse may result in electric shock or fire. To prevent the equipment from being damaged, the following rules should be observed for installation, use and maintenance. **Read the following safety instructions before operating the instrument.**

POWER

POWER CORD: Use only the power cord specified for this equipment and certified for the country of use. If the power (mains) plug is replaced, follow the wiring connections specified for the country of use. When installing or removing the power plug, **hold the plug, not the cord.**

The power cord provided is equipped with a **3-prong grounded plug (a plug with a third grounding pin)**. This is both a safety feature to avoid electrical shock and a requirement for correct equipment operation. If the outlet to be used does not accommodate the 3-prong plug, either change the outlet or use a grounding adapter.

FUSES: Replace fuses only with those having the required current rating, voltage and specified type such as normal blow, time delay, etc. **DO NOT** use makeshift fuses or short the fuse holder. This could cause a shock or fire hazard or severely damage the instrument.

OPERATION

CAUTION

DO NOT OPERATE WITH COVERS OR PANELS REMOVED. Voltages inside the equipment consist of line (mains) that can be anywhere from 100-240VAC.

DO NOT OPERATE WITH SUSPECTED EQUIPMENT FAILURES. If any odor or smoke becomes apparent turn off the equipment and unplug it immediately. Failure to do so may result in electrical shock, fire or permanent damage to the equipment. Contact the factory for further instructions.

DO NOT OPERATE IN AN EXPLOSIVE ATMOSPHERE: Operating the equipment in the presence of flammable gases or fumes **constitutes a definite safety hazard**. For equipment designed to operate in such environments the proper safety devices must be used such as dry air or inert gas purge, intrinsic safe barriers and/or explosion-proof enclosures.

DO NOT IMPEDE THE CHAMBER FROM VENTING EXCESS PRESSURE. The humidification and dehumidification systems are open loop systems that pump external air into the chamber. If the chamber is not allowed to vent, pressure could build up and cause serious damage to the chamber.

IF YOUR UNIT INCLUDES OPTIONAL LIQUID NITROGEN COOLING CAPABILITIES, REVIEW ALL SAFETY INFORMATION IN THE LIQUID NITROGEN SAFETY ADDENDUM.

DO NOT USE IN ANY MANNER NOT SPECIFIED OR APPROVED BY THE MANUFACTURER: Unapproved use may result in damage to the equipment or present an electrical shock or fire hazard.



I. INFORMATIONS DE SÉCURITÉ IMPORTANTES

CONSIGNES DE SÉCURITÉ

L'équipement décrit dans ce manuel est conçu et fabriqué pour fonctionner dans les limites de conception définies. Toute mauvaise utilisation peut entraîner un choc électrique, des dommages à l'équipement ou un incendie. Pour éviter ces événements indésirables, les règles suivantes doivent être respectées pour l'installation, l'utilisation et la maintenance. Lisez les consignes de sécurité suivantes avant d'utiliser l'instrument. Conservez ces instructions dans un endroit sûr pour référence future.

PUISSANCE

CORDON D'ALIMENTATION : Utilisez uniquement le cordon d'alimentation spécifié pour cet équipement et certifié pour le pays d'utilisation. Si la prise d'alimentation (secteur) est remplacée, suivez les connexions de câblage spécifiées pour le pays d'utilisation. **Lors de l'installation ou du retrait de la fiche d'alimentation, tenez la fiche et non le cordon.**

Le cordon d'alimentation fourni est équipé d'une fiche de mise **à la terre à 3 broches standard des États-Unis (une fiche avec une troisième broche de mise à la terre)**. Il s'agit à la fois d'une fonction de sécurité pour éviter les chocs électriques et d'une exigence pour le bon fonctionnement de l'équipement. Ne pas désactiver ou retirer la broche de mise à la terre. Si la prise à utiliser n'est pas compatible avec la fiche à 3 broches, changez la prise ou utilisez un adaptateur de mise à la terre.

FUSIBLES: Remplacez les fusibles uniquement par des fusibles ayant le courant nominal, la tension et le type spécifiés tels que fusion normale, temporisation, etc. **N'UTILISEZ PAS** de fusibles de fortune ou ne court-circuitez pas le porte-fusible. Cela pourrait entraîner un risque d'électrocution ou d'incendie ou endommager gravement l'instrument.

OPÉRATION

MISE EN GARDE

NE PAS UTILISER AVEC LES COUVERCLES OU LES PANNEAUX RETIRÉS. Les tensions à l'intérieur de l'équipement consistent en une ligne (secteur) pouvant aller de 100 à 240 VCA.

NE PAS UTILISER AVEC DES PANNES D'ÉQUIPEMENT SOUPÇONNÉES. Si une odeur ou de la fumée apparaît, éteignez l'équipement et débranchez-le immédiatement. Le non-respect de cette consigne peut entraîner un choc électrique, un incendie ou des dommages permanents à l'équipement. Contactez l'usine pour plus d'instructions.

NE PAS OPÉRER DANS UNE ATMOSPHÈRE EXPLOSIVE: Utilisation de l'équipement en présence de gaz, de poudres ou de fumées inflammables **CONSTITUE UN DANGER DÉFINITIF POUR LA SÉCURITÉ**. Pour les équipements conçus pour fonctionner dans de tels environnements, les dispositifs de sécurité appropriés doivent être utilisés, tels qu'une purge d'air sec ou de gaz inerte, des barrières à sécurité intrinsèque et/ou des boîtiers antidéflagrants.

NE PAS EMPÊCHER LA CHAMBRE D'ÉVACUER L'EXCÈS DE PRESSION. Les systèmes d'humidification et de déshumidification sont des systèmes à boucle ouverte qui pompent l'air extérieur dans la chambre. Si la chambre n'est pas autorisée à se ventiler, la pression peut s'accumuler et causer de graves dommages à la chambre.

SI VOTRE UNITÉ COMPREND DES CAPACITÉS DE REFROIDISSEMENT À L'AZOTE LIQUIDE EN OPTION, VÉRIFIEZ TOUTES LES INFORMATIONS DE SÉCURITÉ DANS L'ADDENDA DE SÉCURITÉ À L'AZOTE LIQUIDE.

NE PAS UTILISER D'UNE MANIÈRE NON SPÉCIFIÉE OU APPROUVÉE PAR LE FABRICANT: Une utilisation non approuvée peut endommager l'équipement ou présenter un risque d'électrocution ou d'incendie.

II. Description of Contents

ETS 5532 and ETS 5533 Environmental chambers. Standard items:

Item	Qty.	Description
Chamber	1	Controlled environment: Acrylic construction, clear and white panels, with doors, glove ports, or iris ports as requested. On the right-hand side is an equipment enclosure which houses user controls and systems for humidity control of the chamber. On the left-hand side, the controlled section of the chamber contains heating and cooling equipment, a circulation fan, and lighting, as ordered. Controller & Sensor Mounted on the front panel of the equipment enclosure, the Controller contains user controls for setting temperature and humidity setpoints and selecting or disabling operating systems. The Controller monitors chamber condition using the calibrated temperature/humidity sensor which enters the controlled section of the chamber through the right-hand wall. Operating systems (if included) A. Heating. Heating elements are mounted inside a panel at the rear of the controlled area. B. Cooling. A thermoelectric cooler is mounted at the rear of the controlled area. C. Humidification. An ETS Model 5482 humidification system is built inside the equipment enclosure and supplies humidified air to the chamber through a 1" tube. D. Dehumidification. Several choices are available to provide dehumidified air and supply to the chamber through ¼" tubing.
5 gallon water tank	1	Plastic tank with valve. NOTE: Shipped with valve inside the tank – remove valve and screw onto tank.
Desiccant column	1	1 Liter column filled with MS desiccant with blue indicating dots.
¼" OD opaque tubing	1	One 10-foot length of tubing is provided to be used for connecting the dehumidification system to the chamber.
Sealing Putty	1	A roll of Mortite for sealing small open ports.
Power Cord	1	AC line cord for connecting the control unit to power.
Pair of Iris Ports	1	If requested, attached to chamber and ready for use.
Pair of Glove Ports	1	If requested, Gloves attach to removable ports. Ready to exchange with Iris ports.
Pair of Glove hands	1	If requested, Extra pair of glove hands for replacement.

III. Setup Guide

Part 1: Unpacking Chamber

The standard Model 5532/5533 Controlled Environment Chamber is shipped in 2 separate cartons. The chamber is double boxed in a double-wall custom carton, and the accessories, shown in Figure 3-1, are packed in a second double-walled carton.

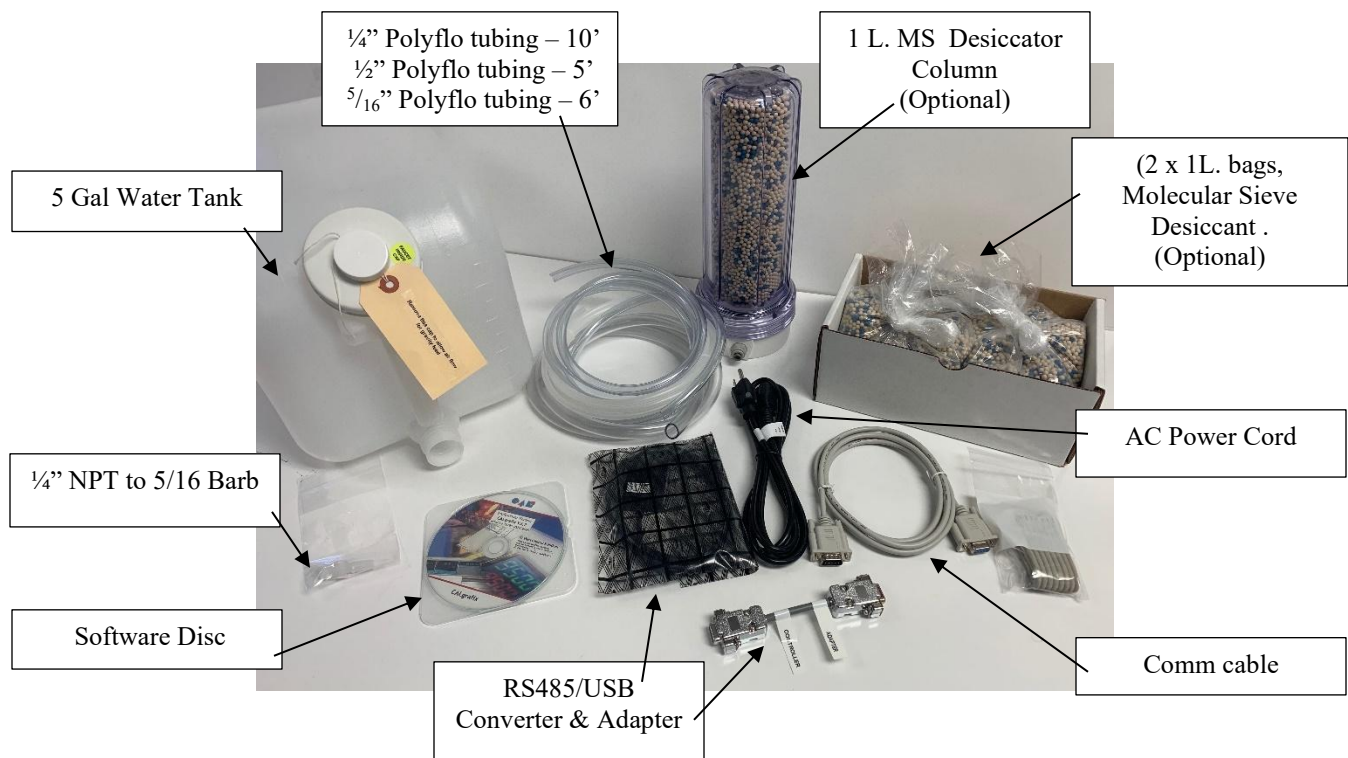


Figure 3-1 Accessory box contents

Chambers with glove ports will also include a spare pair of gloves.

Upon receipt, immediately inspect the cartons for any visible damage. If any shipping damage is noticed, unpack the chamber and inspect it for damage. Take pictures of any abnormalities observed.

Save all cartons for safe transport of the system in the future. If the Chamber has to be returned to ETS for any reason, the original packaging must be used.

NOTE: Report any damage immediately to the common carrier delivering the System and to ETS. All damage claims must originate from the recipient. Failure to report damage in a timely manner may result in the claim not being paid. ETS will not be responsible for damaged or lost components if not reported to ETS within 30 days of shipment.

Part 2: Connect all Operating Systems


Tools Needed: Adjustable wrench

IMPORTANT!

REVIEW ALL SAFETY PRECAUTIONS
AND INSTRUCTIONS

Switches OFF. AC unplugged.

Verify that the POWER switch and all controller switches on the front panel are in the OFF (0) position.

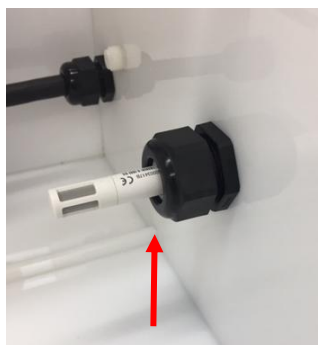

Step 1 – Verify sensor position

The sensor should be inserted into the environmental chamber from the equipment enclosure until the metal portion extends beyond the fitting as shown.

If necessary, loosen the nut, adjust the sensor position, and then tighten the nut to hold the sensor in place.



Valve Closed



Valve Open

Step 2 – Set Up Humidifier Water source

To maximize the life of the transducer, use an appropriate water supply (deionized water or distilled water).

If using the water tank, position the valve as shown, verify open/closed operation.



Step 3. – Position Tank

Select a convenient location for the water tank, above the chamber surface level.



Step 4. – Attach Water Source to Chamber

Use the ¼" clear tubing to route water to the "Water In" connection on the chamber.

Make sure both ends of the tubing are securely inserted into the connectors.


Open the water tank valve.



DESICCANT DEHUMIDIFICATION SYSTEM (if included)

Step 5 – Set Up Dehumidifier

A clear plastic column of Molecular Sieve ("MS") sits outside the chamber. The desiccant column may be placed on top of, or next to, the Chamber.

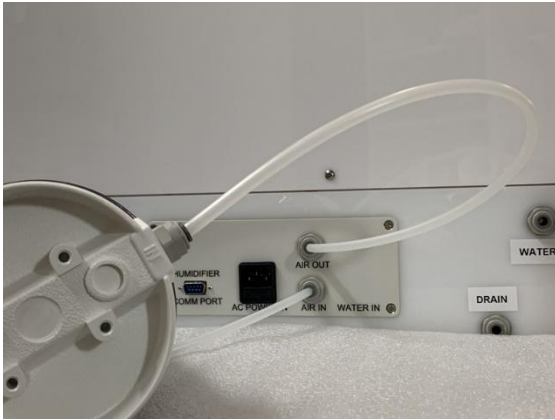

Step 6 – Connect air tubes

Use ¼" O.D. tubing to connect the desiccant column to the chamber through the quick-connect fittings on the right side of the Chamber.

Connect Column IN to chamber AIR OUT.

(Column IN/OUT labels are on the bottom of column.)



Connect Column OUT to chamber AIR IN.

After connecting both tubes, return the column to upright position and set it alongside the chamber.




SELF-REGENERATING DEHUMIDIFICATION SYSTEM (if included)

Step 7 – Connect Air In supply

The (1) self-regenerating dehumidification system is and (2) overpressure monitoring system are installed internally. Overpressure monitoring will cut power to the dehumidification system if an over-pressure event is detected.

Connect the facility compressed air supply using ¼" tubing (included) to the **Air In** connection. Firmly insert the ¼" Polyethylene tubing into the fitting. The **Air Out** connection is not utilized for this system.

The air pressure from the Compressed Air / Dry Gas source must be in the range of **50 – 100 psi (90 psi recommended)** for safe and effective operation.

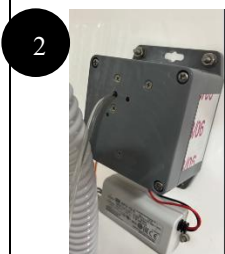



Step 8 – Prepare to power up

Make sure this power switch is OFF.



Step 9 – Connect main power

1. Plug in the AC cable to the chamber.
2. Plug in the to an electrical outlet.

Note: Using 10A 120VAC 60Hz service.

(AC Power IN will be clearly labeled if a different voltage is in use.)

IV. Quick Start Guide

IMPORTANT!

REVIEW ALL SAFETY PRECAUTIONS
AND INSTRUCTIONS

Step 1 – Switches OFF

Before applying power to the system, verify that all controller switches on the front panel are in the OFF (0) position.




Step 2 – Turn on Power

Turn on Power (I = ON). Do NOT turn on “COOL” switch on front panel.



Step 3 – Set Temperature and RH

Set your set points by pressing and holding the * key and using increase/decrease buttons to adjust to your desired target condition.

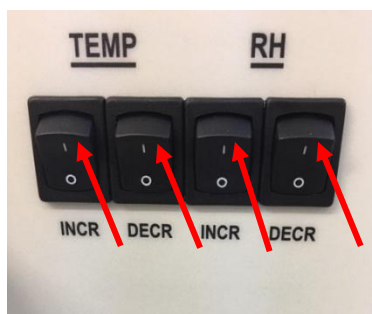

Step 4 – Turn on Operating Systems

Flip ONLY the operating system needed to the ON position (I = on):

Heating: TEMP INCR ON
Cooling: TEMP INCR & DECR ON
Humid.: RH INCR ON
Dehumid.: RH DECR ON

V. Operation Guide

The Model 5532/5533 Microprocessor Controlled Environmental Chamber is a completely integrated system providing a controlled environment for testing, fabricating or storage.

The controllers, pumps and the humidifier operating system are housed in a separate Equipment Enclosure on the right side of the Chamber as shown below.

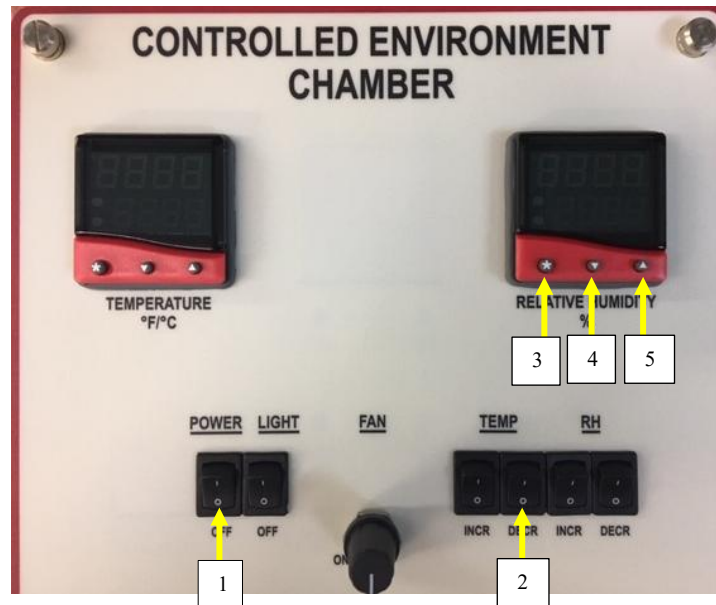


Figure 5-2: Model 5532 Controlled Environment Chamber

All systems are accessible via the removable white acrylic panel on the right side of the Chamber. Controllers and operating systems that are available with the Model 5532/5533 are as follows:

1. Microprocessor Temperature Controller (multi ramp/soak, set point display; Std.)
2. Model 556 Temperature compensated RH and Temperature Sensor (Std.)
3. Integrated 500 Watt Heater (Std.)
4. Choice of Cooling Systems: (5473-150W - Std.)
 - a. Model 5473-(150W-400W) 510-1360 BTU Thermoelectric Cooling System
 - b. Model 5463 Liquid CO₂ Cooling System
 - c. Model 5466 Liquid Nitrogen Cooling System
5. Microprocessor Humidity Controller (multi ramp/soak, set point display, Std)
6. Model 5482 Ultrasonic Humidification System (Std)
7. Choice of Dehumidification Systems:
 - a. Model 5471 Desiccant/Pump Dehumidification System, 2.5 lbs. of MS Desiccant (Std.)
 - b. Model 5478 Self-Regenerating Dehumidification System (requires external air compressor or house air at 50-100 psi)
 - c. Model 5465 Dry Gas Dehumidification System
8. CALGrafix Computer Software/Interface Package. Allows remote monitoring, charting and reprogramming of the Microprocessor Controllers from a PC (Std).

Controller



Button	Description	Functionality
1	POWER	This switch disconnects all power going to the Chamber Systems. "I" is "ON", "O" is "OFF".
2	COOL (TEMP DECR)	Allows the user to manually disable the Temperature Cooling System. Pushing this switch to the ('0') is "Off"
3, 4, 5	TEMPERATURE CONTROLLER	The temperature controller includes a * button along with down and up arrows to allow selecting and changing settings.

Operating Systems: Heating

The Model 5532/5533 Chamber contains two (2), 250 Watt Electric Heaters (500 Watts total). The Heater is an integral part of the Chamber, located inside the chamber workspace. Since the Heater is already installed, there is no additional set-up required.

After turning "on" the "TEMP INCREASE" switch, the Microprocessor Temperature Controller governs the operation of the Heating System.

1. The Microprocessor will determine the amount of heat needed to maintain the desired set point. Power will be applied to the heater as a series of time proportioned pulses. .
2. The electric heating element and indicator light are located on white metal panel containing the Fans located on the inside, rear wall (the heater is hidden from view).
3. The Heater Indicator Light will illuminate to indicate when power is being applied to the Heater. (The Microprocessor LEDs will illuminate at the same time.)

A thermal safety switch (also on the panel) is connected in series with the heater. Power to the heater will be cut off if the temperature within the chamber exceeds 135°F (58°C)

Operating Systems: Cooling

The Model 5532/5533 Chamber includes the *ETS Model 5473-150W, 510 BTU/hr Thermoelectric Cooling System* as standard. Other cooling options include:

ETS Model 5473-300	1020 BTU/hr Thermoelectric Cooling System
ETS Model 5463	Liquid Carbon Dioxide Cooling System
ETS Model 5466	Liquid Nitrogen Cooling System

The Model 5473 Series of thermoelectric cooling systems utilize the Peltier effect to reduce the temperature of a large heat sink. One or more 110 cfm fans circulate the air within the workspace through the heat sink to continuously reduce the temperature. Very precise temperatures can be maintained with this system by the microprocessor temperature controller. The standard Model 5473-150 is capable of removing up to 510 BTU/hr (150 Watts/hr) from the Chamber. The Model 5473-300 is capable of removing up to 1020 BTU/hr (300 Watts/hr) from the Chamber.

Operating Systems: Humidification

The Model 5532/5533 accomplishes humidification using an ETS Model 5482 Ultrasonic Humidification System. The Humidifier is located inside the Chamber Equipment Enclosure. To access the Humidifier, remove the large white acrylic panel on the right hand side of the Chamber. The panel is held in place with eight (8) Phillips screws.

The Ultrasonic Humidifier is an open loop system that produces a fine water mist through ultrasonic action. A blower fan forces the fine mist from the humidifier into the chamber.

Operating Systems: Dehumidification

The Model 5532/5533 Chamber includes a Desiccant-Pump Dehumidification System as standard. As an option, the ETS Model 5478 Self-Regenerating Dehumidification System is available.

Desiccant/Pump Dehumidification System

The Dehumidification System is closed loop designed to reduce the relative humidity in the Model 5532/5533 Chamber to less than 10%. When paired with the Microprocessor Humidity Controller, the humidity inside the chamber can be controlled, without disturbance, to within +/- 0.2% RH of the set point at the sensor.

The Dehumidification System includes a small air pump (located inside the Chamber Control Cavity), 1 liter of Molecular Sieve ("MS") desiccant in a clear plastic column (sits outside the chamber), and ¼" O.D. tubing to interface the drying column with the pump. The tubing

connects to the chamber through the quick-connect fittings on the right side of the Chamber. The desiccant column may be placed on top of, or next to, the Chamber.

The air pump is already connected to the Model 5532/5533 Microprocessor Humidity Controller. All the user needs to do is connect the Desiccant Column into the system using the provided ¼" OD Tubing.

The desiccant removes any moisture that is in the air. This dried air is then forced back into the chamber.

VI. Calibration and Maintenance

Calibration

The only components of the Model 5532/5533 Controlled Environment System that require periodic calibration are the Microprocessor Controllers and Temperature/Humidity sensor. It is recommended that this section be returned to ETS for calibration. However, system calibration can be checked by the user.

Place the sensor of a known, calibrated temperature and humidity meter as close as possible to the System sensor. After the temperature and humidity within the chamber has stabilized, compare the readings. They should be within the combined tolerance of the two measuring devices.

To return equipment to ETS for calibration or repair it is first necessary to obtain a RMA number. **Call 215-887-2196.**


To remove the Control Panel and sensor, first turn off and unplug the chamber, then gain access to the operating systems compartment by removing the 8 screws holding the right side panel in place. Unplug all connections to the front Control Panel. Loosen the 4 captive screws holding the Control Panel and remove. Remove the sensor by loosening the compression fitting inside the controlled environment section. Return both the Control Panel and sensor to ETS.

Maintenance

The Model 5532/5533 requires very little maintenance. The procedures in this manual should be followed to ensure trouble free operation.

Humidifier Preventive Maintenance – Regular Maintenance - Every 500 hours of operation.

Perform regular preventive maintenance (flush out the basin) after every 500 hours of operational use. Any time the humidifier is not being used for more than a week, drain the water completely out of the basin.

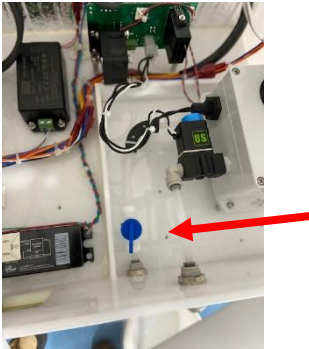POWER OFF!	Step 1 – Turn off power! Unplug the AC line cord from the AC outlet of controller to remove power from the humidifier.
	Step 2 – Prepare and drain the basin The water will need to be drained out of the basin first. First, place the unit above a small container at the end of the drain valve to collect drained water. Turn the drain valve to ON position [tap turned parallel to the valve body]. Wait for the water to stop flowing before proceeding.




Step 3 – Turn on power

Empty the small container, and place back at the end of the drain; leave the drain valve open. Turn the chamber ON, set RH setpoint to 100%, switch on the RH increase to apply constant power to the humidifier.



Step 4 – Flush out contaminants from the basin.

Turn the 5-gallon water tank valve to the “ON” position. Let fresh water from the tank flush out contaminants in the basin through the drain valve for 2 minutes.

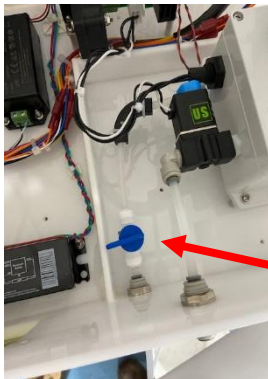

Step 5 – Return the System to operation

Turn the drain valve to the “OFF” position [tap turned to 90° from the valve body]. The humidifier is ready to be placed back into operation again.

Humidifier Maintenance – Transducer Maintenance, If needed.

The humidification system utilizes an ultrasonic transducer that has an estimated lifespan of at least 3000 hours of run time after which it may need to be replaced. The ultrasonic transducer is not user replaceable. Please contact ETS for assistance with repair/replacement.

Dehumidification Maintenance – As needed

In a desiccant system, the dehumidification desiccant blue indicator dots turn a darker blue to indicate that the molecular sieve material needs to be changed.

VII. Troubleshooting.

Troubleshooting Guide		
Problem	Possible Cause	Corrective Action
Chamber not operating – controllers not illuminated.	Is the main chamber fuse blown?	Check the main fuse near the AC cord: If blown, replace with the same type fuse
	Some other power problem?	Check AC Power Source - Is it switched OFF? Is a breaker tripped?
Insufficient Humidification – unable to reach target RH% level:	Water basin is empty.	Verify water tank latch is turned to “on” position
		Verify basin drain is turned to “closed” position
		Verify water tank has at least 1 gallon of water
		Verify that water tank is located above basin
		Verify that air relief on water tank is open/loose
	Preventative maintenance not on schedule.	Perform maintenance on the humidifier unit. Do this regularly.
	Dew point limits have been reached	Is the temperature high or low, creating a dew point limit?
	Too much backpressure causing the humidification system to not operate properly	Verify enough chamber/glove box venting.
Unstable RH% level		Review controller setup
Insufficient Dehumidification – unable to reach target dry level	Problem with Desiccant System	Are all tubing connections properly connected? Is desiccant used up? Is there a moisture source being used in the chamber?
	Problem with Dry gas dehumidification	Are all tubing connections properly connected? Is the tank empty or gas source turned off? Is there proper pressure and flow?

	Problem with self-regenerating dehumidification:	Are all tubing connections properly connected? Is a compressor connected to the system? Is there proper pressure and flow from the compressor? Is the compressor spraying water due to lack of maintenance?
Insufficient Heating	Low ambient temperature or leaks into chamber. Heater panel not activating when needed.	Review controller setup and connections. Is the lab ambient temperature lower than expected? Are there leaks from outside? Does the light on the heater panel illuminate when the temperature microprocessor green indicator light illuminates?
Insufficient Cooling	Heat source in the chamber or leaks from ambient. External air flow insufficient to remove heat.	Review controller setup and connections. Is the lab ambient temperature higher than expected? Is there a heat load inside the chamber? Is the fan on the back outside of the chamber running? Is there 2' of clearance behind the chamber with good air flow?
Excessive noise / vibration	Related to motors or fans.	Check for sources such as an item contacting fan blades.

VIII. System Specifications

Chamber:

Construction:

Enclosure: 3/8" (6mm) clear & white acrylic

Polished welded museum quality seams

Doors: 1/2" (13mm) clear acrylic, 1/2-turn latches

Front door: 32 x 14", Side door: 4 x 14"

Seal: 1/4" (6mm) Closed-cell, non-setting gasket

Access:

Ports: 2x 1/4" (6mm) Quick disconnect (right side)

1x 1" (25.4mm) Hose barb (right side)

1x 1 1/2" (31.4mm) cable pass through (left side)

Gloves: (optional) 2x 8" (20.3cm) ports.

Circulating Fans: 2x 110 cfm (2 x 3115 L/min),
ON/VARIABLE/OFF control

Lighting: Enclosed LED lighting

Dimensions:

5532: External: 54"Wx25.5"Dx22"H (137x65x56cm)

Internal usable working space: 39"Wx20"Dx21"H
(99x51x54 cm)

5533: External: 54"Wx31.25"Dx22"H (137x80x56cm)

Internal usable working space: 39"Wx26"Dx21"H
(99x66x53 cm)

Weight: 5532: 120 lbs. (54 kg) 5533: 136 lbs.(62 kg)

Operating Range:

Humidity: <10 to >98%

Temp: (Std) 64-135°F (18-55°C)

Controllers:

Control unit: 2x CAL 9500, multiple ramp/soak cycle
(Optional): 1x CAL 9500, 3rd parameter, linear input

Display Resolution: 0.1% / 0.1°

Recorder Out: 0-1vdc Analog.

Computer Interface: MODBUS proto. via RS485-
PC with W98, NT, W Pro 7, or higher

Software: (optional) CALGrafix PC chart & data collect.

Sensor:

Probe, Temperature Compensated RH/Temp

RH: Capacitive film, 0 to 100%, non-condensing

Accuracy: $\pm 1.5\%$ RH @ 72°F (22 °C)

Temperature: RTD, -40 to +140 °F (-40 to 60 °C)

Accuracy $\pm 0.5\%$ °F ($\pm 0.3\%$ °C)

Operating Systems: (If included)

Dehumidify: Desiccant/Pump system, including
2x 1 L bag of Molecular Sieve desiccant .

Humidify: Ultrasonic micro-mist generation

Capacity: DI water line or 5 gal (19 L) ext. tank

Flow rate – 0.06 cfm (max)

Heat: Electric, 500W

Cool: 150W Thermoelectric (Standard), 300W (Optional)

Power:

Voltage: 115 VAC, 10 Amps, 50/60Hz

Optional: 230 VAC, 5 Amps, 50/60 Hz

Dual North American GFI AC outlet inside chamber

NOTE: The entire humidity range cannot be obtained at all temperatures due primarily to dew point considerations as shown below.

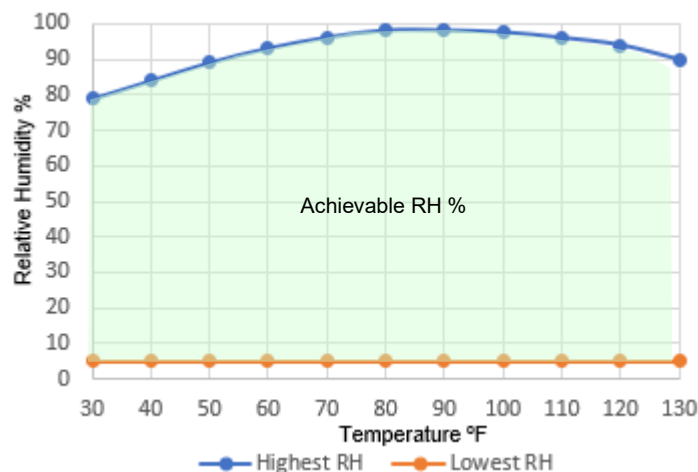


Figure 5-1: Model 5532/5533 Dew point performance chart

Controller and operating system Options that are available with the Model 5532/5533 are as follows:

- 1) Microprocessor Temperature Controller (multi ramp/soak, set point display; Std.)
- 2) Model 556 Temperature compensated RH and Temperature Sensor (Std.)
- 3) Integrated 500 Watt Heater (Std.)
- 4) Choice of Cooling Systems:
 - a) Model 5473-(150W) 510BTU Thermoelectric Cooling System (Std)
 - b) Model 5473-(400W) 1360 BTU Thermoelectric Cooling System
 - c) Model 5463 Liquid CO₂ Cooling System
 - d) Model 5466 Liquid Nitrogen Cooling System
- 5) Microprocessor Humidity Controller (multi ramp/soak, set point display, Std)
- 6) Model 5482 Ultrasonic Humidification System (Std)
- 7) Choice of Dehumidification Systems:
 - a) Model 5471 Desiccant/Pump Dehumidification System, 1 L. of MS Desiccant (Std.)
 - b) Model 5478 Self-Regenerating Dehumidification System (requires external air compressor or house air at 50-100 psi)
 - c) Model 5465 Dry Gas Dehumidification System
- 8) CALGrafix Computer Software/Interface Package. Allows remote monitoring, charting and reprogramming of the Microprocessor Controllers from a PC (Std).

IX. Warranty

Limited Warranties. Seller warrants that all goods manufactured and delivered hereunder shall (a) conform to any samples, drawings, specifications or other written documents provided to Seller by Buyer, or approved by Buyer to Seller and (b) be free from all defects in workmanship and material. Buyer's sole remedy against Seller for breach of either of the specifically mentioned warranty shall be the repair or replacement, at Seller's sole option, of the defective workmanship or material. Seller expressly disclaims all other warranties, express and/or implied, including but not limited to those of merchantability and fitness for a particular purpose.

In no event shall Seller be liable, under either warranty or otherwise, to Buyer in excess of the purchase price of the products paid to Seller by Buyer. In no event shall Seller be liable for any loss or damage arising directly or indirectly from the use of the product or for consequential or incidental damages. Seller's specified warranties will expire and lapse (i) for renewable items (such as gloves, iris ports and desiccants), sixty (60) days from date of shipment and (ii) for all standard equipment and otherwise nonrenewable items, one year from date of shipment.